




multiSUB miniRAPIDE

The **miniRAPIDE** makes imaging your gels simple. With its UV transparent base there is no need to handle the gel directly, perfect for teaching or space restricted labs.

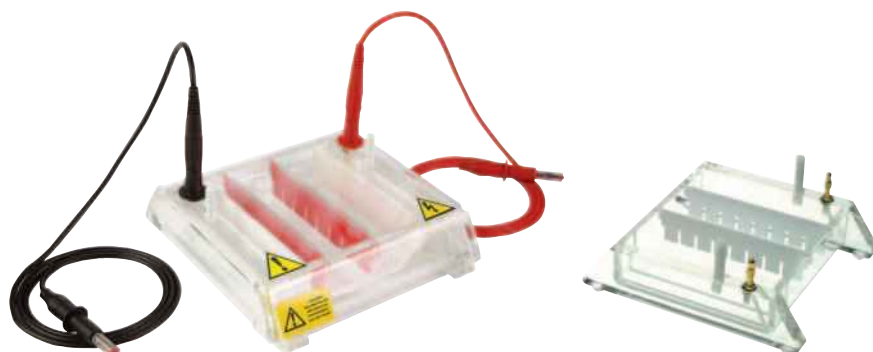

KEY FEATURES

miniRAPIDE is intended for agarose gel electrophoresis

- UV Transparent for direct gel imaging
- Low buffer volumes for low cost running
- Multichannel pipette compatible combs for speed

TECHNICAL SPECIFICATIONS

Gel dimensions,	10 x 8cm (W x L)
Unit dimensions	15 x 15 x 4cm (W x D x H)
Max. sample capacity	40 samples
Buffer volume	50ml
Combs available : No. of samples Thicknesses	1, 4, 8, 12, 16, 20 1, 1.5mm

ORDERING INFORMATION

FMMS10 miniRapid, 10 x 8cm, 2 x 8 sample combs 1.5mm and casting dams

FMMS-DAM miniRapid Casting dams, pk/2
RPW0.2 Replacement Platinum Wire 0.2mm – 50cm

Combs



FMMS-1-1, (330µl)



FMMS-4-1, (90µl)



FMMS-8-1, (40µl)



FMMS-12-1, (25µl)



FMMS-16-1, (15µl)



FMMS-20MC-1, (10µl)

For 1.5mm combs, replace -1 with -1.5 in the ordering code e.g. FMMS-4-1.5.

multiSUB-4

multiSUB-4 is a compact system capable of running over 1200 samples simultaneously by stacking up to 4 horizontal gels.

Each multiSUB-4 is supplied with 4 gel trays and 8 combs as standard. Two double-sided comb and three tray length formats, 8x6, 8x12 and 8x18cm are also available. These multichannel-compatible combs and gel plate configurations are compatible with microplates and thermal cycler blocks to ensure rapid loading of restriction digests and PCR products by 8-channel pipette.



KEY FEATURES

- Separates a maximum 1200 samples in as little as 15 minutes in 4 stacked gel trays
- Double-sided 1.5mm thick combs allow more sample volume to be loaded into each well
- Three gel tray options available in 8x6, 8x12 and 8x18cm (WxL) sizes for maximum flexibility
- Optional Flexicaster

TECHNICAL SPECIFICATIONS

Gel dimensions (w x l)	8 x 6cm, 8 x 12cm, 8 x 18cm
Unit dimensions	11 x 35 x 16cm
Max. sample capacity per 18cm tray	with 1cm run length: 306 with 2cm run length: 144 with 3cm run length: 72
Buffer volume	200, 400, 600 or 800ml (for 1, 2, 3 or 4 gel trays resp.)
Combs available: No. of samples Thicknesses	1, 8, 12, 18 DuoCombs 1, 1.5mm

ORDERING INFORMATION

CSL-MULTISUB4 multiSUB-4 multi-level Gel Chamber, includes 4x 18cm UV Trays, 8x 18/8 Sample 1.5mm Combs (Tape UV Trays to seal)

CSL-MULTISUB4EXCAS Multisub-4, as above but with External Caster for 4 gels

MSUB4-12/1-1 12/1 Sample 1mm Combs for multiSUB-4

MSUB4UV6 multiSUB-4 tray 8 x 6cm

MSUB4-18/8-1 18/8 Sample 1mm Combs for multiSUB-4

MSUB4UV12 multiSUB-4 tray 8 x 12cm

MSUB4-12/1-1.5 12/1 Sample 1.5mm Combs for multiSUB-4

MSUB4UV18 multiSUB-4 tray 8 x 18cm

MSUB4-18/8-1.5 18/8 Sample 1.5mm Combs for multiSUB-4